

SUPINE TO PRONE IN ONE MOTION







Zing MPS, Supine, MPS Tilt Table(TT), Supine Tilt Table(TT), Vertical, & Prone

The Zing line of multi and single position standers were designed to fit the widest range of children and accommodate positioning challenges few others can. The Zing line features independent leg abduction up to 30 degrees on each model. The units can be configured simply, with padded strapping system to ease transfers and aid multi-user environments. But when required, individual positioning and support components can be added.



The / Ing

True Multi-Position Stander

- The MPS base allows three positions without changing pads, flipping foot plates or moving the tray
- Provides a smooth transition from a flat-to-load transfer to supine, vertical and 20 degrees prone
- Frame will only go supine and upright with Prone Position Lockout engaged

True Hip Pivot Leg Abduction

- · Pivots at the hips for proper leg alignment in hips, knees and feet
- Feet maintain angle with the rest of the legs no wedges needed
- Up to 30 degrees of abduction per leg
- Promotes musculoskeletal development for infants, toddlers and kids

Plenty of Growth

- Size 1 fits children from infancy to 44"
- · Can easily grow with the child and adapt to their changing needs
- 70 pound weight capacity

Multi-Use Trav

- The MPS, Supine, MPS TT and Supine TT feature a multi-use tray that is angle adjustable and swings away
- Tray can remain attached during transfers
- More activities are possible by swinging tray away during standing therapy

Advantage

Accommodating Contractures

- with contractures
- or both sides

Hygienic Covers

- Easy to wash
- Keeps stander from getting soiled

Unique Lift Mechanisms

- hand control

Remarkably Modular

- complex needs
- as needed

• No special modifications are needed to accommodate children

• Using the multi-adjustable foot plates, multi-adjustable knee pads, and the independent calf pads, it is possible to adjust for contractures in one

• Available for the pelvic, upper body, lateral, head and tray areas

· Especially beneficial if stander is used during feeding

• Dual controls are available allowing angle adjustment by foot pedal or

 Gas Spring Lockout prevents an inadvertent change of position, especially important when around other children

· Can be configured as a simple stander or fit a child with the most

• Designed to grow with the child – options can be added

· Great for a multi-user facility or school





Zing MPS, Supine, MPS Tilt Table(TT), Supine Tilt Table(TT), Vertical, & Prone

SUPINE TO PRONE



123














SPACE RACE







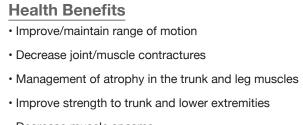

- Decrease muscle spasms
- Promote proper muscle skeleton development
- Improve/maintain bone integrity/skeletal development
- Lessen/manage the progression of scoliosis
- Manage pressure (ulcers) through changing positions
- Improve bowel function and regularity
- Aid in kidney and bladder functions
- Strengthen cardiovascular system and build endurance
- Improve circulation
- Reduce swelling in lower extremities

Zing Prone



Hygienic Cover Patterns





SEA CREATURES



FANCY PANTS ELEPHANTS






ALL SMILES ANIMALS



FEATHERS IN FLIGHT

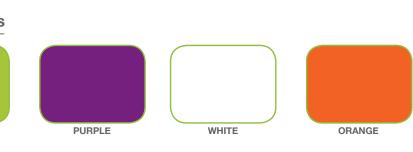





Zing MPS, Supine, MPS Tilt Table(TT), Supine Tilt Table(TT), Vertical, & Prone

SPECIFICATIONS


additional specifications at easystand.com

The Zing MPS is the only multi-position stander that can seamlessly transition from supine the child out and reposition footThe Zing Supine features a comfortable 32" transfer height and provides all the positioning and support features of the MPS.The Zing MPS TT is identical to the MPS from a positioning standpoint, but features larger, wider support pads and options for paddedThe Zing Supine TT features the larger pads and padded strapping options of the MPS TT and positioning from flat to load supine to vertical.The Zing Vertical is a single position stander that features a padded Quick Strap to ease upright transfers. It's also one of the only vertical standers to offerThe Zing Prone is a single position stander that features a pads and options for padded to vertical.The Zing Supine TT features the larger pads and padded strapping options of the MPS TT and positioning from flat to load supine to vertical.The Zing Vertical is a single position stander that features a padded Quick Strap to ease upright transfers. It's also one of the only vertical standers to offerThe Zing Prone is a single position stander that features a to vertical.The Zing Supine TT features the larger pads and padded strapping options of the MPS TT and positioning from flat to load supine to vertical.The Zing Supine TT features the position stander that features a padded Quick Strap to ease upright transfers. It's also one of the only vertical standers to offerThe Zing Prone is a single transfer the only vertical standers to offerThe Zing Prone is a single transfer the duice standerThe Zing Prone is a single transfer the duice standerThe Zing Prone is a single to vertical.The Zing Prone is a single to positio
plates, pads or trays. transferring and securing children. anatomically correct independent leg abduction. degrees.
Weight Limit 70 lbs. (32kg)
Height Range up to 44" (112cm)
Transfer Height 32" (81 cm)
Pelvic Support to Foot Plate Range 7"-22" (18-56cm) 7"-22" (18-56cm) 7"-22" (18-56cm) 7"-22" (18-56cm)
Upper Body Support Height Range 4" (10cm) 4" (10cm) 4" (10cm) 4" (10cm) 4" (10cm)
Knee Pad Depth Range 5" (13cm) 5" (13cm) 5" (13cm) 5" (13cm) 5" (13cm)
Weight of Basic Unit 55 lbs (25kg)
Frame Footprint 20"x28" (51x71cm)
30° Independent Leg AbductionAvailableAvailableAvailableAvailable
Supine 0-90° 0-90° 0-90° n/a
Prone 0°-20° n/a 0°-90°
HCPCS E0641 E0638 E0643 E0638 E0638
Positions Supine, Vertical, Prone Supine, Vertical, Prone Supine, Vertical, Prone Supine, Vertical, Prone Vertical, Prone
Tray Adjustable, Swing-Away Adjustable, Swing-Away Adjustable, Swing-Away Adjustable, Swing-Away Adjustable, Swing-Away
SupportIndividual Positioning Components, Padded StrappingIndividual Positioning Components, Padded Strapping
Dual ControlAvailableAvailableAvailableAvailable
Gas Spring Lockout Available Available Available Available Available
Transfer Aid Flat-to-Load Flat-to-Load Flat-to-Load Flat-to-Load flat-to-Load flat-to-Load flat-to-Load
Tool Free AdjustmentYesYesYesYes



А

FITTING THE ZING

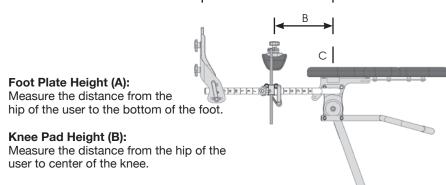

Centering the Hip (C): The user's hip should be centered in line with the pivot point of the stander ensuring proper hip alignment when using abduction, adjust support components up or down based hip position.

ZING MPS (Multi-Position Stander) Options

The MPS (multi-position stander) base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver with Prone Position Lockout



Comes with four swivel locking casters, foot operated gas spring, mast receiver with prone position lockout.

PA5520 Zing MPS Base



PA5530 Mast Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head supports.



PA5532 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required to do leg abduction.

PA5612 Dual Control

Allows both foot and hand control of the gas spring lift. Must order push handle A5594.

PA5614 Gas Spring Lockout

Locks gas spring in any position from lying down to fully upright and prone positions.



PA5608 Directional Locking Caster Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

PA5536 Platform Foot Plate

Adjustable in height. Foot holders and foot straps can be added for more positioning of the feet.



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm),

PA5540, PA5542, PA5544 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots

for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) belt with a D-Ring Velcro closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.





PA5560 Multi-Adjustable Calf Pads 3.5" PA5678 Multi-Adjustable Calf Pads 4.5" Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable up and down with knee pad and fore/aft 4" (10cm) independently. Must order


multi-adjustable knee pads.









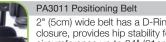

PA5610 Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).

PA3028, PA3029 Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5562 and PA5564.

PA3030, PA3032 Hygienic Covers

Washable slip on cover available in six different patterns. For form to fit pads PA5566 and PA5570.



0

70 lbs



PA5566 Form to Fit Pads-5"Hx10"W PA5570 Form to Fit Pads-5"Hx12"W

Form to fit pad can be custom fitted to the user. Pad sizes: PA5566-5"Hx10" (13x25cm) PA5570-5"Hx12"W (13x30cm). Not recommended with leg abduction mast.

PA5574 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm). For planar pads only.

















PA5562 Planar Pads-5"Hx7"W PA5564 Planar Pads-5"Hx9"W

Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13x18cm), PA5564-5"Hx9"W (13x23cm).

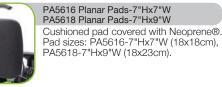

PA5568 Form to Fit Pads-5"Hx10"W PA5620 Form to Fit Pads-7"Hx10"W

Form to fit pad with built in lateral support that can be custom fitted to the user. Pad size for PA5568 is 5"Hx10"W (13x25cm) and PA5620 s 7"Hx10"W (18x25cm).

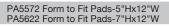

Form to fit pad with built in lateral support that can be custom fitted to the user. Pad size for PA5572 is 5"Hx12"W (13x30cm) and PA5622 is 7"Hx12"W (18x30cm).

PA5576 Lateral Supports

Keeps the upper torso in symmetrical alignment. Adjustable from 7-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar pads only.



PA5580

PA5578 Lateral Supports with Elbow Stop Works as a lateral support and elbow rest Great for the supine position. For planar pads only.



Lateral Supports w/Elbow Stop/Arm Rest

Works as a lateral support, elbow rest and armrest. Great for the supine, upright and prone position. For planar pads only.

PA5630 Elbow Stop with Arm Rest

Great for loading and bringing to upright. Works with both planar and form to fit upper body supports. For planar and form to fit pads.



A PS PT50086 X-Style Chest Vest-9.5"Lx9"W PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.



PA3028, PA3029, PA3045, PA3046 Hygienic Covers



PA5576.

PY5628.

1

PA3037 Hygienic Cover

Washable slip on cover available in six different patterns. For planar pads PA5562, PA5564, PA5616, PA5618.



PA3031, PA3033, PA3047, PA3048 Hygienic Covers Washable slip on cover available in six different patterns. For form to fit pads PA5568, PA5572, PA5620, PA5622.



PA5602 Clear Swing Away Tray-11"Lx19"W Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5604 Forearm Wings-Black Molded Tray

Wings give additional area to rest arms when in standing position.

PA5606 Forearm Wings-Clear Acrylic Tray

Wings give additional area to rest arms

when in standing position.



PY5626 Head Support-5"Hx8"W PY5628 Head Support-6"Hx10"W

PA5586 Form to Fit Headrest

upper body pad is 8" (20cm).

PA3035, PA3036 Hygienic Covers

different patterns. For PY5626 and

Cushion size is 5"Hx8"W (13x20cm) 6"Hx10"W (15x25cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from seat is 24" (61cm)

Cushion size 5"Hx10"W (13x25cm). Depth

range is -2.5" (6cm) to +1.5" (4cm) from

backrest. Maximum height from top of

Washable slip on cover available in six

Washable slip on cover available in six

PA5596 Black Molded Tray-11"Lx19"W

Tray is accessible in all positions and

rom supine to upright and prone

Tray size is 11"Lx19"W (28x48cm).

provides anterior support as user moves

positions. Adjustable in depth and height.

different patterns. For PA5586.



PA3038 Hygienic Padded Cover Washable slip on cover is padded and available in six different patterns.



PNG50042 Hand Grips Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction

cups. Clear tray required.



Makes pushing the stander easier for the

P80809 Tool Pouch




Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.



P80246 Angle Locator

Helps therapists monitor standing progress. Attaches to the mast. Can be used on leg abduction.

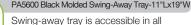

positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

ZING Supine Options

The Supine base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



PA5522 Zing Supine Base Comes with four swivel locking casters, foot operated gas spring and a mast receiver.



PA5530 Mast Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head supports.



PA5532 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required to do leg abduction.

PA5612 Dual Control

Allows both foot and hand control of the gas spring lift. Must order push handle PA5594.

PA5614 Gas Spring Lockout

Locks gas spring in any position from lying down to fully upright and supine positions.



PA5608 Directional Locking Caster Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

PA5536 Platform Foot Plate

Adjustable in height. Foot holders and foot straps can be added for more positioning of the feet.



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi -/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm).

PA5540, PA5542, PA5544 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps

measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) belt with a D-Ring Velcro closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.





PA5560 Multi-Adjustable Calf Pads 3.5" PA5678 Multi-Adjustable Calf Pads 4.5" Great for adapting to knee contractures.

Gives support from the back of the calf. Adjustable up and down with knee pad and fore/aft 4" (10cm) independently. Must order multi-adjustable knee pads.












(5cm) belt. Fits circumference up to 24" (61cm).

PA3028, PA3029 Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5562 and PA5564.

PA3030, PA3032 Hygienic Covers

Washable slip on cover available in six different patterns. For form to fit pads PA5566 and PA5570.

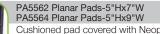



PA5566 Form to Fit Pads-5"Hx10"W PA5570 Form to Fit Pads-5"Hx12"W

Form to fit pad can be custom fitted to the user. Pad sizes: PA5566-5"Hx10" (13x25cm) PA5570-5"Hx12"W (13x30cm). Not recommended with leg abduction mast.

PA5574 Hip Supports


Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm). For planar pads only.

PA3011 Positioning Belt

2" (5cm) wide belt has a D-Ring Velcro closure, provides hip stability for user. Fits circumference up to 24" (61cm).

PA5610 Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2"

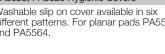









70 lbs



















PA5562 Planar Pads-5"Hx7"W PA5564 Planar Pads-5"Hx9"W

Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13x18cm), PA5564-5"Hx9"W (13x23cm).

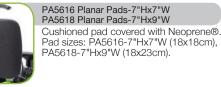

PA5568 Form to Fit Pads-5"Hx10"W PA5620 Form to Fit Pads-7"Hx10"W

Form to fit pad with built in lateral support that can be custom fitted to the user. Pad size for PA5568 is 5"Hx10"W (13x25cm) and PA5620 s 7"Hx10"W (18x25cm).

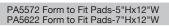

Form to fit pad with built in lateral support that can be custom fitted to the user. Pad size for PA5572 is 5"Hx12"W (13x30cm) and PA5622 is 7"Hx12"W (18x30cm).

PA5576 Lateral Supports

Keeps the upper torso in symmetrical alignment. Adjustable from 7-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar pads only.



PA5578 Lateral Supports with Elbow Stop Works as a lateral support and elbow rest Great for the supine position. For planar pads only.



PA5580

armrest. Great for the supine, upright and prone position. For planar pads only.


PA5630 Elbow Stop with Arm Rest


Great for loading and bringing to upright. Works with both planar and form to fit upper body supports. For planar and form to fit pads.



PT50086 X-Style Chest Vest-9.5"Lx9"W A PS PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.



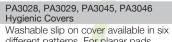

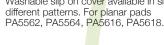



1

PA3031, PA3033, PA3047, PA3048 Hygienic Covers Washable slip on cover available in six different patterns. For form to fit pads PA5568, PA5572, PA5620, PA5622.

PY5626 Head Support-5"Hx8"W

PY5628 Head Support-6"Hx10"W

PA5586 Form to Fit Headrest

upper body pad is 8" (20cm).

PA3035, PA3036 Hygienic Covers

different patterns. For PY5626 and

Washable slip on cover available in six

Washable slip on cover available in six

PA5596 Black Molded Tray-11"Lx19"W

Tray is accessible in all positions and

rom supine to upright and prone

Tray size is 11"Lx19"W (28x48cm).

provides anterior support as user moves

positions. Adjustable in depth and height.

different patterns. For PA5586.

Cushion size is 5"Hx8"W (13x20cm)

6"Hx10"W (15x25cm). Depth range is +1"

Maximum height from seat is 24" (61cm)

Cushion size 5"Hx10"W (13x25cm). Depth

range is -2.5" (6cm) to +1.5" (4cm) from

backrest. Maximum height from top of

to -2.75" (+2 to -7cm) from backrest.

PA5576.

PY5628.

PA3037 Hygienic Cover



PA5602 Clear Swing Away Tray-11"Lx19"W Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5604 Forearm Wings-Black Molded Tray

Wings give additional area to rest arms when in standing position.





PA5606 Forearm Wings-Clear Acrylic Tray Wings give additional area to rest arms when in standing position.



PA3038 Hygienic Padded Cover

Washable slip on cover is padded and available in six different patterns.



PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.



PA5594 Push Handle

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.

P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.



P80246 Angle Locator

Helps therapists monitor standing progress. Attaches to the mast. Can be used on leg abduction.

PA5600 Black Molded Swing-Away Tray-11"Lx19"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

Zing MPS Tilt Table (TT) Options

The MPS (multi-position stander) base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver with Prone Position Lockout



Comes with four swivel locking casters, foot operated gas spring, mast receiver with prone position lockout.

PA5520 Zing MPS TT Base




PA5632 Mast Mast is the base for support and positioning components.



PA5634 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required to do leg abduction.





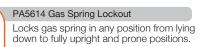

Allows both foot and hand control of the

gas spring lift. Must order push handle

PA5608 Directional Locking Caster

PA5612 Dual Control

PA5594.



Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

PA5536 Platform Foot Plate

Adjustable in height. Foot holders and foot straps can be added for more positioning of the feet.



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi -/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm).

PA5540, PA5542, PA5544 Foot Holders

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm). PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5"

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5644, PA5646, PA5648 Leg Support Planar Pads Cushioned Neoprene® pad for leg support. Pad sizes: PA5644-4"Hx6"W (10x15cm) PA5646-6"Hx6"W (15x15cm) PA5648-8"Hx6"W (20x15cm)

Washable slip on cover available in six different patterns. For planar pads PA5644, PA5646, PA5648.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) belt with a D-Ring Velcro closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.















PA5666, PA5668, PA5670 Knee Strap

4"x6" (10x15cm) pad attached to a 1.5" (4cm) belt. Fits circumference up to 15" (38cm). Adjustable in height. For leg supports PA5644, PA5646, PA5648

PA5556 Multi-Adjustable Knee Pads-4"Hx3.5"W PA5558 Multi-Adjustable Knee Pads-4"Hx4"W

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5556 width 4"-9" (10x23cm) and PA5558 4.5"-9" (11x23cm). With a depth range of 5" (13cm).

PA5560 Multi-Adjustable Calf Pads 3.5" PA5678 Multi-Adjustable Calf Pads 4.5"

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable up and down with knee pad and fore/aft -1"-4" (-2.5-10cm) independently. Must order multi-adjustable knee pads.

PA5638, PA5640, PA5642 Planar Pad

Cushioned pad covered with Neoprene®. Pad sizes: PA5638-14"Hx14"W (36x36cm)

PA5640-16"Hx14"W (41x36cm) PA5642-18"Hx14"W (46x36cm)

PA3063, PA3064, PA3065 Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5638, PA5640, PA5642.

PA5636 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm).

PA3062 Positioning Strap

2" (5cm) wide belt has a D-Ring Velcro closure, provides hip stability for user. Fits circumference up to 24" (61cm).



70 lbs

12



















PA5660, PA5662, PA5664 Hygienic Covers







PA5672 Pelvic Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).



PA5652 Lateral Supports-4"-8"Range PA5654 Lateral Supports-7"-11"Range Keeps the upper torso in symmetrical alignment. PA5652 Adjustable from 4"-8" (10-20cm), PA5654 Adjustable from 7"-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm).




PA5674 Upper Body Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).



PT50086 X-Style Chest Vest-9.5"Lx9"W PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.



PA3039 Hygienic Covers for Lateral Supports Washable slip on cover available in six different patterns. For lateral supports PA5652 and PA5654.

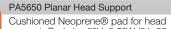

support. Pad size 8"Hx8.5"W (21x22cm).



PY5626 Head Support-5"Hx8"W PY5628 Head Support-6"Hx10"W

Cushion size is 5"Hx8"W (13x20cm), 6"Hx10"W (15x25cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from seat is 24" (61cm).

PA5586 Form to Fit Headrest



Cushion size 5"Hx10"W (13x25cm). Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).

PA5676 Headrest Mounting Bracket

For use when mounting an after market head support.



Washable slip on cov<mark>er a</mark>vailable in six different patterns. For head supports PA5650, PY5626, PY5628, PA5586,




PA5596 Black Molded Tray-11"Lx19"W

Tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height. Tray size is 11"Lx19"W (28x48cm).


PA5600 Black Molded Swing-Away Tray-11"Lx19"W



Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

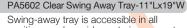


positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

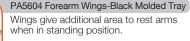

PA5606 Forearm Wings-Clear Acrylic Tray Wings give additional area to rest arms when in standing position.



PNG50042 Hand Grips

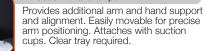



Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.



P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.



P80246 Angle Locator

Helps therapists monitor standing progress. Attaches to the mast. Can be used on leg abduction.

Zing Supine Tilt Table (TT) Options

The Supine base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



Comes with four swivel locking casters, foot operated gas spring, and a mast receiver.

PA5522 Zing Supine TT Base



PA5632 Mast Mast is the base for support and




positioning components.



PA5634 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required to do leg abduction.





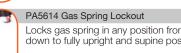

PA5612 Dual Control

PA5594.

Locks gas spring in any position from lying down to fully upright and supine positions.


Allows both foot and hand control of the

gas spring lift. Must order push handle



PA5608 Directional Locking Caster locks one caster while leaving the other three free to swivel.

PA5536 Platform Foot Plate

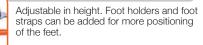



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi -/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm).

PA5540, PA5542, PA5544 Foot Holders

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm). PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5644, PA5646, PA5648 Leg Support Planar Pads Cushioned Neoprene® pad for leg support. Pad sizes: PA5644-4"Hx6"W (10x15cm) PA5646-6"Hx6"W (15x15cm) PA5648-8"Hx6"W (20x15cm)

PA5660, PA5662, PA5664 Hygienic Covers Washable slip on cover available in six different patterns. For planar pads PA5644, PA5646, PA5648.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) belt with a D-Ring Velcro closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.















PA5666, PA5668, PA5670 Knee Strap

4"x6" (10x15cm) pad attached to a 1.5" (4cm) belt. Fits circumference up to 15" (38cm). Adjustable in height. For leg supports PA5644, PA5646, PA5648

70 lbs

12

PA5556 Multi-Adjustable Knee Pads-4"Hx3.5"W PA5558 Multi-Adjustable Knee Pads-4"Hx4"W

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5556 width 4"-9" (10-23cm) and PA5558 4.5"-9" (11-23cm). With a depth range of 5" (13cm).

PA5560, PA5678 Multi-Adjustable Calf Pads

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable up and down with knee pad and fore/aft -1"-4" (-2.5-10cm) independently. Must order multi-adjustable knee pads.

PA5638, PA5640, PA5642 Planar Pad

Cushioned pad covered with Neoprene®. Pad sizes: PA5638-14"Hx14"W (36x36cm)

PA5640-16"Hx14"W (41x36cm) PA5642-18"Hx14"W (46x36cm)

PA3063, PA3064, PA3065 Hygienic Covers

Washable slip on cover available in six different patterns. For planar pads PA5638, PA5640, PA5642.

PA5636 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm).

PA3062 Positioning Strap

2" (5cm) wide belt has a D-Ring Velcro closure, provides hip stability for user. Fits circumference up to 24" (61cm).





















Great for transporting stander. Directionally

PA5672 Pelvic Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).



PA5652 Lateral Supports-4"-8"Range PA5654 Lateral Supports-7"-11"Range Keeps the upper torso in symmetrical

alignment, PA5652 Adjustable from 4"-8" (10-20cm), PA5654 Adjustable from 7"-11 (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm).




PA5674 Upper Body Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).



PT50086 X-Style Chest Vest-9.5"Lx9"W PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.



PA3039 Hygienic Covers for Lateral Supports Washable slip on cover available in six different patterns. For lateral supports PA5652 and PA5654.

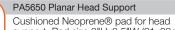

support. Pad size 8"Hx8.5"W (21x22cm).



PY5626 Head Support-5"Hx8"W PY5628 Head Support-6"Hx10"W


Cushion size is 5"Hx8"W (13x20cm), 6"Hx10"W (15x25cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from seat is 24" (61cm).


PA5586 Form to Fit Headrest



Cushion size 5"Hx10"W (13x25cm). Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).

PA5676 Headrest Mounting Bracket

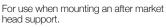

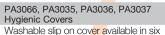

different patterns. For head supports PA5650, PY5626, PY5628, PA5586,




PA5596 Black Molded Tray-11"Lx19"W

Tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height. Tray size is 11"Lx19"W (28x48cm).


PA5600 Black Molded Swing-Away Tray-11"Lx19"W





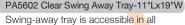

positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

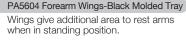



PA5606 Forearm Wings-Clear Acrylic Tray Wings give additional area to rest arms when in standing position.




PA3038 Hygienic Padded Tray Cover

Washable slip on cover for tray is padded and available in six different patterns.



PNG50042 Hand Grips

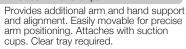

PA5594 Push Handle



P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.





P80246 Angle Locator

Helps therapists monitor standing progress. Attaches to the mast. Can be used on leg abduction.



Zing Vertical Options

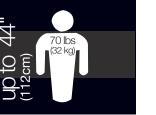

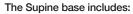

- Four Swivel Locking Casters
- Mast Receiver



(25x10cm).



Widi P.

Zing Prone Options

The Supine base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



Comes with four swivel locking casters, foot operated gas spring, and a mast receiver.

PA5522 Zing Prone Base



PA5530 Mast Mast is the base for support and positioning components.



PA5532 Mast with Leg Abduction

Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required to do leg abduction.



PA5725 Dual Control

Allows both foot and hand control of the gas spring lift. Must order push handle PA5720.



PA5614 Gas Spring Lift Lockout Locks gas spring in any position from lying

down to fully upright and supine positions.



PA5608 Directional Locking Caster Great for transporting stander. Directional locks one caster while leaving the other three free to swivel.

PA5680 Platform Foot Plate



Adjustable in height. Foot holders and foot 🕅 straps can be added for more positioning of the feet.



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm),

PA5540, PA5542, PA5544 Foot Holders

Jse with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm). PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L Velcro® and D-ring adjustment holds feet

in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) belt with a D-Ring Velcro closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.


PA5688 Multi-Adjustable Knee Pads-4"Hx3.5"W PA5690 Multi-Adjustable Knee Pads-4"Hx4"W

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5688 width 4"-9" (10-23cm) and PA5690 4.5"-9" (11-23cm). With a depth range of " (13cm).

PA5726 Knee Strap

1"x6" (10x15cm) pad attached to a 1.5" (4cm) belt. Fits circumference up to 15" (38cm). Adjustable in height. For PA5688 and PA5690.



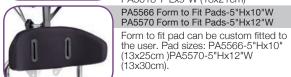







PA3011 Positioning Strap 2" (5cm) wide belt has a D-Ring Velcro

(9x9cm). For planar pads only.

PA5574 Hip Supports

closure, provides hip stability for user. Fits circumference up to 24" (61cm).

Keeps the hips in symmetrical alignment.

PA5610 Padded Positioning Strap


5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).

PA3092 Quick Strap

2" (5cm) wide belt has a airline style buckle for guick closure, provides hip stability for user. Fits circumference up to 24" (61cm).

PA3093 Padded Quick Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt has a airline style buckle for quick closure. Fits circumference up to 24" $(61 \, \text{cm})$

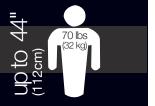































Altimate Medical, Inc. 262 West First Street Morton, MN 56270 USA Tel: 800.342.8968 • 507.697.6393 Fax: 877.342.8968 • 507.697.6900 info@easystand.com • easystand.com

Standing technology should only be used under the guidance of a physician with recommendations for standing program protocol and any medical precautions. Standing programs should be monitored by the attending therapist. AMI maintains a policy of continual product improvement and reserves the right to change features, specifications, and prices without prior notification. Check with AMI for latest information. FORM ZINGCAT101716 REV A Copyright © 2016 Altimate Medical, Inc. All rights reserved. Printed in the USA.